

EnSURE™ Touch

Sanitation Verification.
Designed for You.



hygiena.com/ensuretouch-fb



EnSURE™ Touch

Today's food safety professionals face significant challenges: complex supply chains, increased audit and inspection scrutiny, and the ever-present risk of contamination and recall. Professionals need to ensure processes are free of contamination and efficiently perform product safety interventions.

EnSURE™ Touch is the next-generation quality monitoring system. It provides rapid and accurate sanitation verification to support your FSMA program. Featuring a 5-inch touch screen, wireless sync technology, and cloud-based software, the EnSURE Touch is designed to adapt to your workplace and provides the data you need for complex multi-location operations, audit and risk management, and recall prevention.

EnSURE Touch Can Help You

- Verify sanitation of surfaces, liquids, and product testing in as little as 15 seconds
- Test and compare the effectiveness of sanitation products and procedures
- Train and monitor sanitation crew and QC/QA employees
- Support SSOP and objectively prove cleaning standards have been met
- Quickly identify when replaceable parts like belts and brushes should be replaced
- Identify problem areas in your facility that need more attention
- Provide on-the-spot training and feedback

In Order To:

- Meet FSMA, GFSI, and ISO standards
- Protect your brand reputation
- Avoid financially damaging lawsuits
- Reduce risk of food poisoning, facility closure, and product recalls
- Lower facility testing costs by bringing testing in-house
- Provide evidence of due diligence - HACCP





Innovation at Your Fingertips

Quick test: Run a test from any plan or location in as few as two touches.

Locations: Create testing points with customized fields and pass/fail limits on the go without needing desktop software.

Plans: Easily group testing locations into plans and randomize test locations to ensure sampling coverage.

Users: Manage user access to differentiate personnel from managers. Enhanced administrator control.

Results & Reports: Easily view testing results and trend analysis.

Search: Quickly find specific testing plans and locations on the fly, limiting time between tests.

Sync: Wirelessly synchronize data to SureTrend™ Cloud.

Calibration: Factory-calibrated and preset checks to stay GMP/ISO-compliant.

Retest: Swipe through failed test locations and instantly run re-tests to show effectiveness of corrective action.

MicroSnap™: Test and show results for organisms like Coliform, *E. coli*, Enterobacteriaceae, and total viable count.*

Enzymes: Test and show results for enzymes like Alkaline Phosphatase and Acid Phosphatase.*

Settings: Customize settings for language, custom fields, test types, and simulate competitor RLU limits.

Specifications: Size: 7" x 3" x 1" • Weight: 0.7 lbs • Memory: 2GB internal + cloud storage • Connectivity: USB or Wi-Fi*
Language: Supports multiple languages • Battery: Rechargeable Lithium-Ion battery with USB-C charging *





Dashboard



Results



Reports



Instrument



The data you need to master risk management all in one place.

SureTrend™ Cloud is an easy-to-use data analysis software that comes free with every EnSURE Touch. Available in both cloud-based and desktop formats, SureTrend Cloud enables users to monitor, track, and trend testing results across multiple facilities, and schedule automatic reports, making risk management easier than ever.


Features & Benefits:

- Communicate department performance and achievement to your team and to management
- Fully interactive dashboard allows for quick data viewing, showing you everything at a glance
- Easily configure one or hundreds of instruments from a single SureTrend Account
- Monitor testing across one or multiple facilities
- Schedule automatic reports delivered directly to your email
- Full online support and resource center only one click away
- Track re-tests to discover opportunities for retraining or identify equipment that may need replacing
- Quickly create charts and graphs of test results with many user configurable reports
- Free upgrades with optional wireless software updates

suretrend.hygiena.com

Features & Benefits

- Advanced photodiode technology and unique chemistry for maximum sensitivity and durability
- Intuitive smart navigation makes it easy for anyone to use
- Functions like a smartphone
- USB or Wi-Fi sync capability
- Supports multiple languages

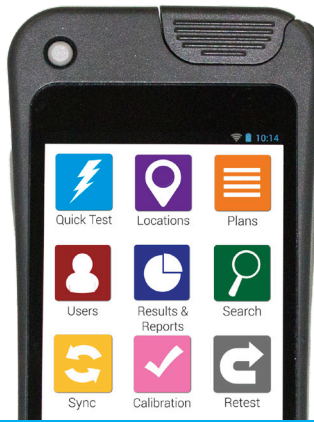

- Responsive 5-inch shatter-proof touch screen
- Train remote teams with built-in screen sharing capability
- Built-in tilt sensor for more accurate readings
- Auto shut-off for smart battery management
- Simulate competitor upper and lower RLU limits

Compatible Products

UltraSnap™ | Surface ATP Test



Catalog No. US2020 | Qty: 100

[Instructional Video Online](#)

SuperSnap™ | High-Sensitivity Surface ATP Test



Catalog No. SUS3000 | Qty: 100

Catalog No. SUS3000X | Qty: 100

[Instructional Video Online](#)

AquaSnap™ TOTAL | Water ATP Test



Catalog No. AQ-100X | Qty: 100

[Instructional Video Online](#)

AquaSnap™ FREE | Water ATP Test



Catalog No. AQ-100FX | Qty: 100

[Instructional Video Online](#)

CalCheck™ | LED Calibration Verification Device



Catalog No. CAL | Qty: 1

[Instructional Video Online](#)

MicroSnap™ | Indicator Organism Tests



Catalog No. MS2-EB | Qty: 100

Catalog No. MS2-ECOLI | Qty: 100

Catalog No. MS2-TOTAL | Qty: 100

Catalog No. MS2-COLIFORM | Qty: 100

[Instructional Video Online](#)

ZymoSnap™ ALP | ALP Pasteurization Verification Test



Catalog No. ZS-ALP-100 | Qty: 100

[Instructional Video Online](#)

CROSS CHECK™ | ACP Cooking Efficiency Test




Catalog No. CX-3000 | Qty: 100

Product Name	Catalog No.	Quantity
EnSURE™ Touch	ETOUCH	1
UltraSnap™ Surface ATP Test	US2020	100
SuperSnap™ High-Sensitivity Surface ATP Test	SUS3000	100
SuperSnap™ High-Sensitivity Surface ATP Test	SUS3000X	100
AquaSnap™ Total Water ATP Test	AQ-100X	100
AquaSnap™ Free Water ATP Test	AQ-100FX	100
MicroSnap™ Coliform & E. coli Enrichment Device	MS1-CEC	100
MicroSnap™ Enhanced Nutrient Broth	MS1-N-BROTH-9ML	100
MicroSnap™ Coliform Detection Device	MS2-COLIFORM	100
MicroSnap™ E. coli Detection Device	MS2-ECOLI	100
MicroSnap™ Enterobacteriaceae Enrichment Device	MS1-EB	100

Product Name	Catalog No.	Quantity
MicroSnap™ Enterobacteriaceae Detection Device	MS2-EB	100
MicroSnap™ Total Enrichment Device	MS1-TOTAL	100
MicroSnap™ Total Detection Device	MS2-TOTAL	100
ZymoSnap ALP Alkaline Phosphatase Test	ZS-ALP-100	100
CrossCheck Acid Phosphatase Test	CX-3000	100
CalCheck LED Calibration Verification Device	CAL	1
ATP Positive Control Kit	CK25	25
ZymoSnap ALP Positive Control Kit	ZS-ALP-PC	5
Small Format Incubator	INCUBATOR	1
Lab Format Incubator	INCUBATOR2	1
BAX System Q7	ASY2018	1
BAX System X5	MCH2007	1

EnSURE Touch is factory calibrated and runs a calibration check on each startup. Recommended calibration control kits available.

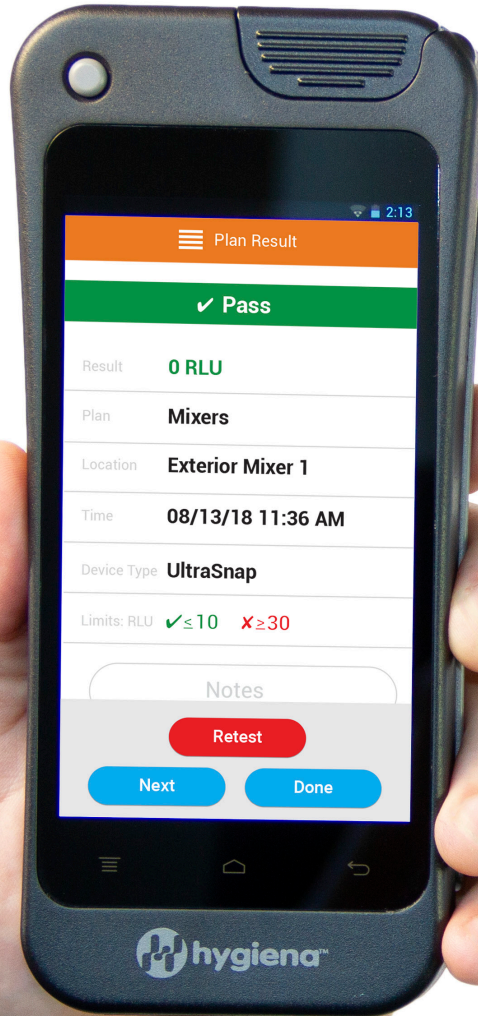

Americas:

Hygiena Headquarters

941 Avenida Acaso
Camarillo, CA 93012
Phone: 1-805-388-8007
Fax: 1-805-388-5531
info@hygiena.com

Hygiena Canada

2650 Meadowvale Blvd Unit 14
Mississauga, Ontario L5N 6M5
Phone: 1-833-494-4362 (Toll-free)
or 1-416-686-7962

Hygiena Mexico, S.A. de C.V.

Calle 3 Anegas 409 Bodega 5, Col. Nueva
Industrial Vallejo, Delegación Gustavo A.
Madero, C.P. 07700, CDMX, México.
Phone: +52-5281-4108 or +52-5281-4146

International:

Hygiena International

8, Woodshots Meadow
Watford, Hertfordshire
WD18 8YU, UK
Phone: +44 (0)1923-818821
Fax: +44 (0)1923-818825

Hygiena (Shanghai) Trading Co., Ltd.

Rm.7K, No.518, Shangcheng Rd.
Pudong New District
Shanghai, China 200120
Phone: +86 21-5132-1081, +86 21-5132-1077,
+86 21-5132-1078
Fax: +86 21-5132-1080

Hygiena Diagnóstica España S.L.

P. I. Parque Plata, Calle Cañada Real 31-35,
41900, Camas, Sevilla, Spain
Phone: +34 954-08-1276

help.hygiena.com

 [YouTube.com/HygienaTV](https://www.youtube.com/HygienaTV)

 [Vimeo.com/HygienaTV](https://vimeo.com/HygienaTV)

